

## IGP 96G

Thermal Cycler Gradient

## IGP 96GR

Thermal Cycler Standard


Thermal cycler is essential laboratory instrument in the field of molecular biology research such as sequencing, gene cloning, gene expression, mutagenesis applications and applied market like drug discovery, agriculture, food industry, etc.



Excellent temperature uniformity




Gradient function

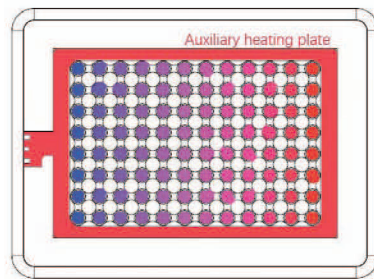

### Features

- High performance long life Peltier and independent control circuits for different heating segments implement accurate temperature control
- Auxiliary heating mechanism diminishes the "edge effect" and enhance the temperature uniformity
- Wide touchdown PCR temperature range (-9.9°C~+9.9°C) and Long PCR time range (-9min59s~+9min59s)
- Gradient temperature setting supported, optimizing temperature easily in the single run
- User friendly interface on 7 inch color touch screen enable you edit the program very simple
- The PCR touch screen pen improves the operation experience and reduces the cross contamination risk

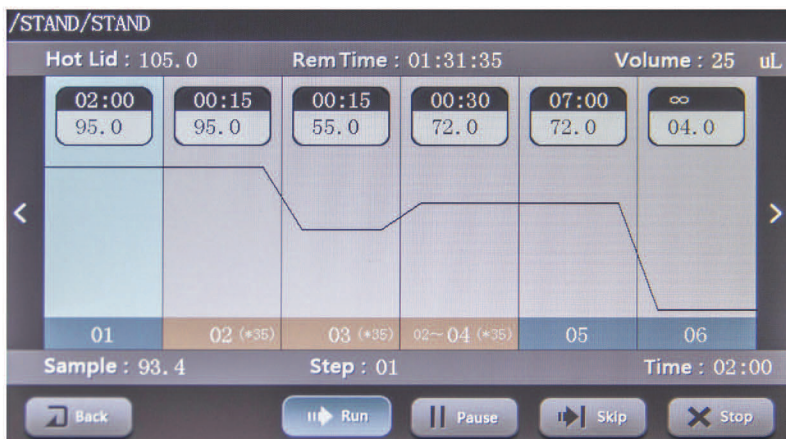

### More features


- Wide options of consumables common PCR tubes, 8-well PCR strips and 96-well PCR plates can be used
- File customization, multi-file storage
- Power-off protection function, automatic program recovery
- Hot lid auto-off function: If the module temperature is lower than 30°C, the hot lid function will automatically turn off



LCD touch panel



Excellent temperature accuracy

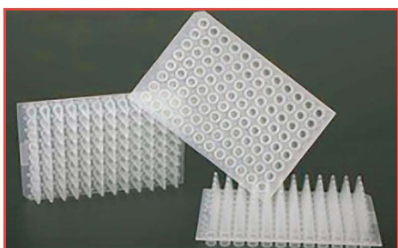

### Specifications

	IGP 96G	IGP 96GR
Sample Capacity	96X0.2mL PCR tube, 8X12 PCR plate or 96 well plate	96X0.2mL PCR tube, 8X12 PCR plate or 96 well plate
Heating Temperature Range /°C	4~105	4~105
Lid Temperature Range /°C	30~110	30~110
Temperature Display Accuracy /°C	±0.1	±0.1
Temperature Display Accuracy@55°C	±0.3	±0.3
Temperature uniformity@55°C	<0.3	<0.3
Max. Heating/Cooling Rate	5°C/Sec	5°C/Sec
Gradient Temperature Setting Range /°C	30~99°C	-
Gradient Range /°C	1~42°C	-
Adapter block material	Aluminum	Aluminum
Display	7" LCD 800x480	7" LCD 800x480
Input	Touch panel	Touch panel
User defined file system	Max. 30 segments 99 cycles max. 16 folder and 16 files each folder	Max. 30 segments 99 cycle max. 16 folder and 16 files each folder
Power off protection	Yes	Yes
Power Supply	100~120V/200~240V,50/60Hz	100~120V/200~240V,50/60Hz
Dimension[WxDxH]	280x370x250 mm	280x370x250 mm
Weight	11kg	11kg